


Energy dependence of the γ -absorption coefficient

(Item No.: P2524215)

Curricular Relevance

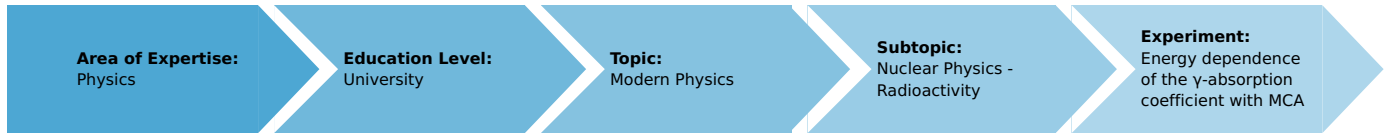

Difficulty



Difficult

Preparation Time



1 Hour

Execution Time



2 Hours

Recommended Group Size



2 Students

Additional Requirements:

- PC

Experiment Variations:

Keywords:

Compton scattering, photoelectric effect, pair production, absorption coefficient, radioactive decay, gamma-spectroscopy

Overview

Short description

Principle

The intensity of γ -radiation decreases when it passes through solid matter. The attenuation can be the result of Compton scattering, the photoelectric effect or the pair production. An absorption coefficient can be attributed to each of the three phenomena. These absorption coefficients, as well as the total absorption, are highly energy-dependent. The energy dependence of the total absorption coefficient for aluminium is measured in the range below 1.3 MeV.



Fig. 1: Experimental set-up.

Equipment

Position No.	Material	Order No.	Quantity
1	Multichannel analyser	13727-99	1
2	measure Software multi channel analyser	14452-61	1
3	Radioactive source Am-241, 370 kBq	09090-11	1
4	Radioactive source Na-22, 74 kBq	09047-52	1
5	Gamma detector	09101-00	1
6	High Precision Power Supply 1.5 kV DC	09107-99	1
7	High-voltage connecting cable	09101-10	1
8	Base plate for radioactivity	09200-00	1
9	Plate holder on fixing magnet	09203-00	1
10	Lab jack, 160 x 130 mm	02074-00	1
11	Vernier calliper stainless steel 0-160 mm, 1/20	03010-00	1
12	Source holder on fixing magnet	09202-00	1
13	Absorption material, aluminium	09029-03	1
14	Screened cable, BNC, l = 750 mm	07542-11	1

Tasks

- Record spectra of the isotopes ^{22}Na and ^{241}Am and determine useful energy windows for the absorption experiment.
- Record the number of incidents during a fixed time in a chosen energy window (integration measurement) in dependence on absorber layer thickness between source and detector for different energy windows.

Setup and Procedure

Set up the experiment as shown in Fig. 1. Before turning on the operating unit for the scintillation counter, connect the high voltage cable correctly to operating unit and photomultiplier and read the instructions in the manual of the gamma-detector. Set the voltage of the operating unit to 1.0 kV. Connect the MCA to the computer's USB port and start the "measure" program. Select the Gauge "Multi Channel Analyzer" and you will receive the start window, as shown in Fig. 2. For recording of spectra you may bring one of the sources in direct vicinity to the detector - for the 370 kBq source you better keep some distance since the counting rate should not exceed 800 cts/s by far.

1. Start with the ^{22}Na source. Select "Spectra recording" (see Fig. 2) and click on "Continue". Set the "Gain" to "Level 1", choose "Channel number" as x-Data (see Fig. 3) and start data recording. The 1275 keV peak should be observable in full at the high energy end of the spectrum, as shown in Fig.3. If not, alter the high voltage setting on the detector operating unit. After adjustment of the high voltage leave it unchanged throughout the experiment.

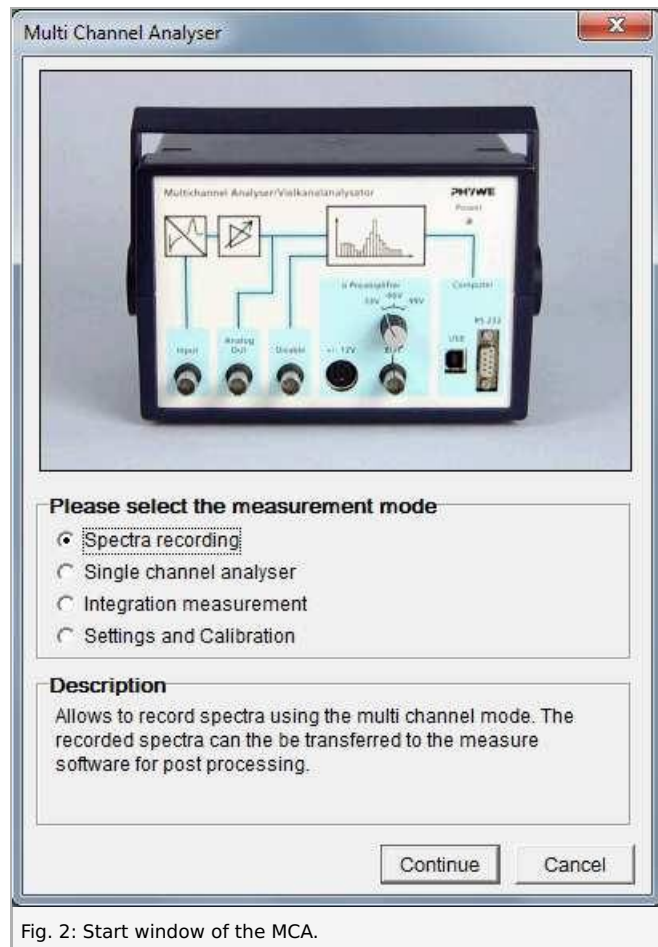
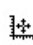


Fig. 2: Start window of the MCA.

2. Note down the channel number ranges of the 511 keV annihilation peak and the 1275 keV γ -peak - the windows you are to determine may be as broad as the photo peaks are. Use the "Mark" tool  on the lower left of the window for channel evaluation. Then set the amplification to "Level 4" and record a spectrum of ^{241}Am - the γ -peak at 60 keV is also suitable for the measurement, as shown in Fig. 4. The chosen window in Fig.4 might be channel 430 to channel 760.

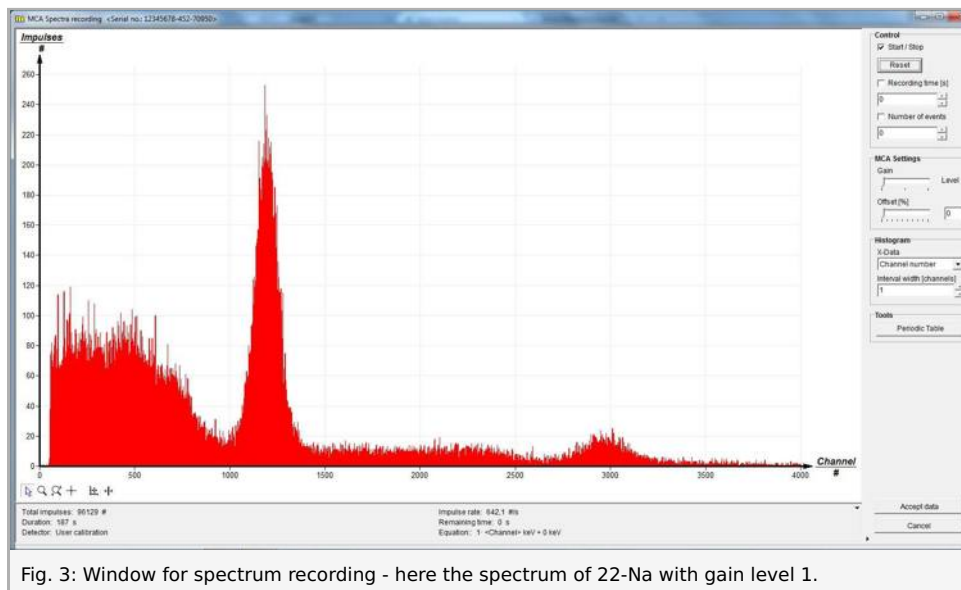


Fig. 3: Window for spectrum recording - here the spectrum of 22-Na with gain level 1.

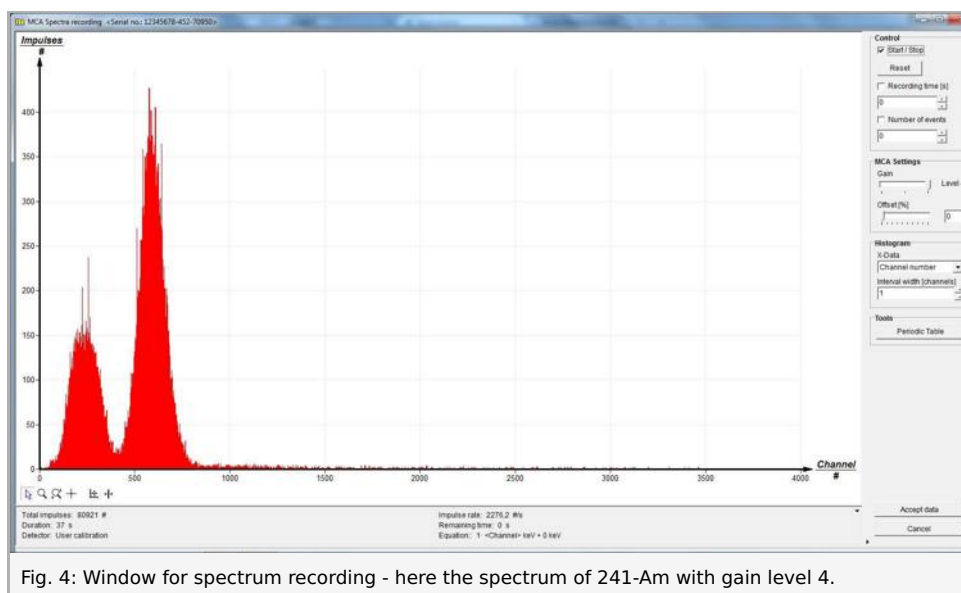


Fig. 4: Window for spectrum recording - here the spectrum of 241-Am with gain level 4.

- Now start the Gauge "Multi Channel Analyzer" again and choose the program part "Integration measurement" from the start window of the MCA (see Fig. 2) and you will receive the integration measurement window, as shown in Fig. 5. Select a source and set the lower and upper limit of the channel to the values you have determined for the peak of interest before. Select a gate time of at least 60 s. Check radiation and photomultiplier dark count rate background in the chosen window to be negligible by performing one measurement without a source and delete the data. Then place the source about at least 15 cm from the detector and make sure, that source and detector are well aligned.

Place aluminium absorber layers in steps of e.g. 5 mm midway between source and middle point of the detector crystal without moving source or detector (placing thin absorbers into the plate holder and thick absorbers onto the supporting blocks) and perform measurements for layers up to 25 mm. The geometry is of importance because scattering on the metal crystal lattice occurs and the counting rate may be altered by the presence of metal even without absorption or Compton-scattering. Change the x-data in the data table of the integration measurement window from number to the appropriate values in mm. Then click on the "Accept data" button. An example of the measurement is shown in Fig. 6 below.

In the measure program use the "Display options" button to set the scaling of the measurement on the "Channels" chart to "logarithmic". A linear dependence of the logarithm of counts on absorber layer thickness should be visible. Save the measurement file to your hard disk drive. Do so for the photo peaks of 1275 keV, 511 keV and 60 keV. (Recording spectra for a defined time, "Accept data" and integrating over the peak with "Mark" tool and "Show integral" button and background subtraction has to be done here for correct intensity values for each absorption layer thickness - the Compton background is the result of higher energetic photons in the detector that get less attenuated.)

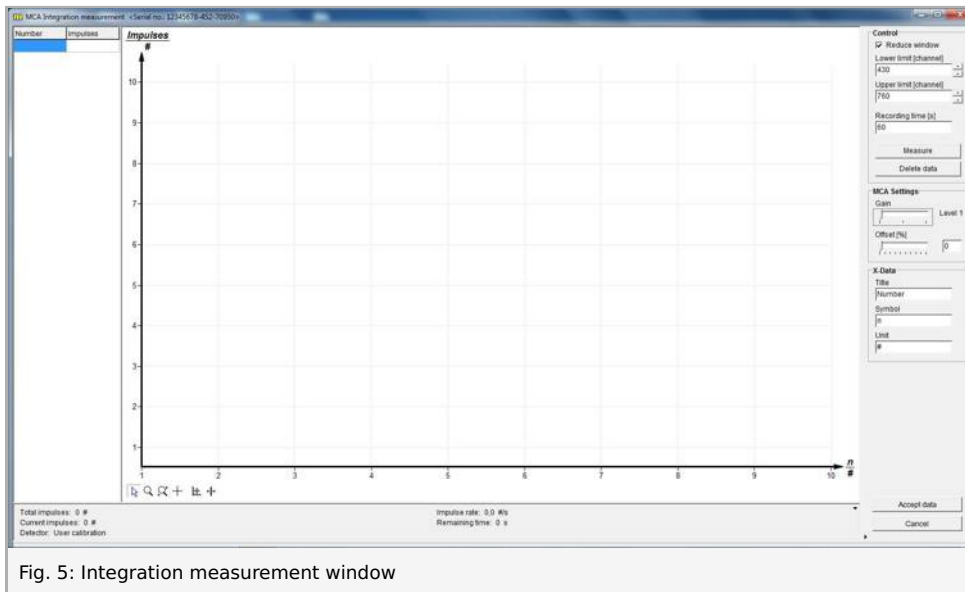


Fig. 5: Integration measurement window

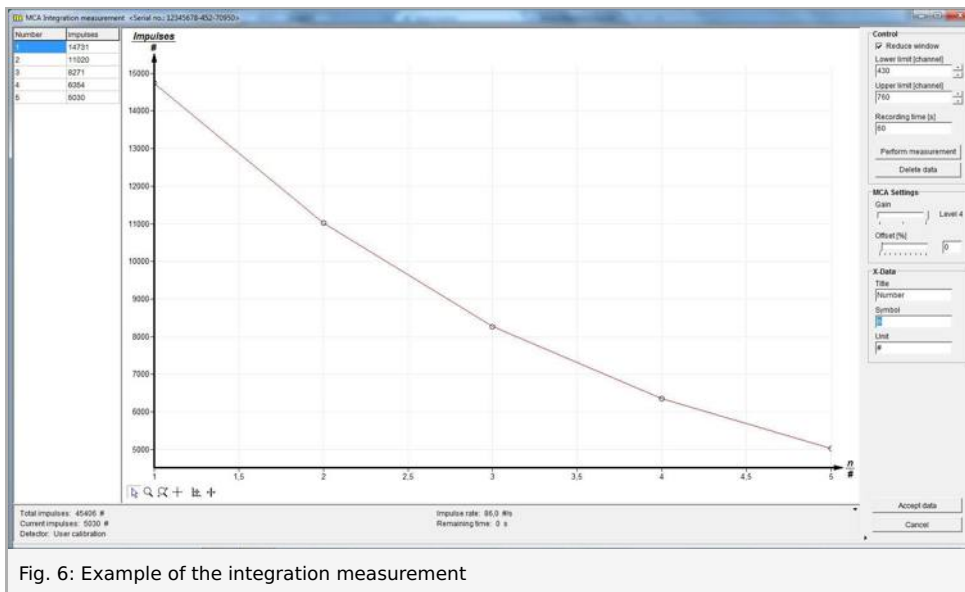


Fig. 6: Example of the integration measurement

Theory and evaluation

Fig. 5 shows the decay schemes of the used nuclides. The proportions of the energy scale are not displayed correctly and the term scheme of ^{237}Np is strongly simplified – in the experiment only the 59.5 keV radiation of ^{237}Np is of importance. It is to be kept in mind that (in the case of ^{137}Ba) the excited states of the daughter nuclides can also disintegrate by inner conversion leading to a strong 32 keV X-ray line. The α -particles, the electrons and positrons are stopped in the sources, the latter giving rise to 511 keV annihilation radiation. Only the generated photons reach the detector.

The γ - and X-ray quanta coming from the sources (α - and β -particles cannot leave covered sources) passing through matter may undergo mainly these reactions:

- They may be diffracted on a crystal lattice without energy loss, which is here therefore not counted for as absorption.
- They may be absorbed by an electron transferring all the energy and momentum to the electron which is called photoelectric effect. In case of condensed matter, thereafter the electron loses its energy in short time ($< 10^{-11}$ s) on a short path (< 1 mm) to the surrounding matter producing mostly heat but maybe also some Cherenkov radiation and a path of ions or excited states of nearby lattice impurities which afterwards recombine partly emitting visible light. The photoelectric effect is dominant for small γ -energies but the cross-section for photoelectric effect strongly decreases for high energies – the likelihood of an electron to swallow a photon completely is only high if the photon has far less energy than $m_e c^2$.
- Else the γ -quanta may interact with electrons exchanging momentum and energy with an electron in an elastic scattering process which is called Compton scattering. In this process, they lose energy and so this process is counted as absorption – the γ -quanta are no longer found in the energy window where they were emitted. This process contributes the most to the total absorption coefficient in the energy range where the photon's energy is comparable to the electron's resting energy $m_e c^2$ (billiard is played with balls of comparable weights).
- If the γ -quanta have more energy than $2 \cdot m_e c^2 = 1022$ keV, with surrounding matter as third partner for momentum conservation, pair production of an electron-positron pair is possible. This process gets more likely with increasing energy – the excess energy separating the reaction products. From the following annihilation of the positron with another electron originate two γ -quanta of 511 keV (or with far less probability three γ -quanta with in average 341 keV each).

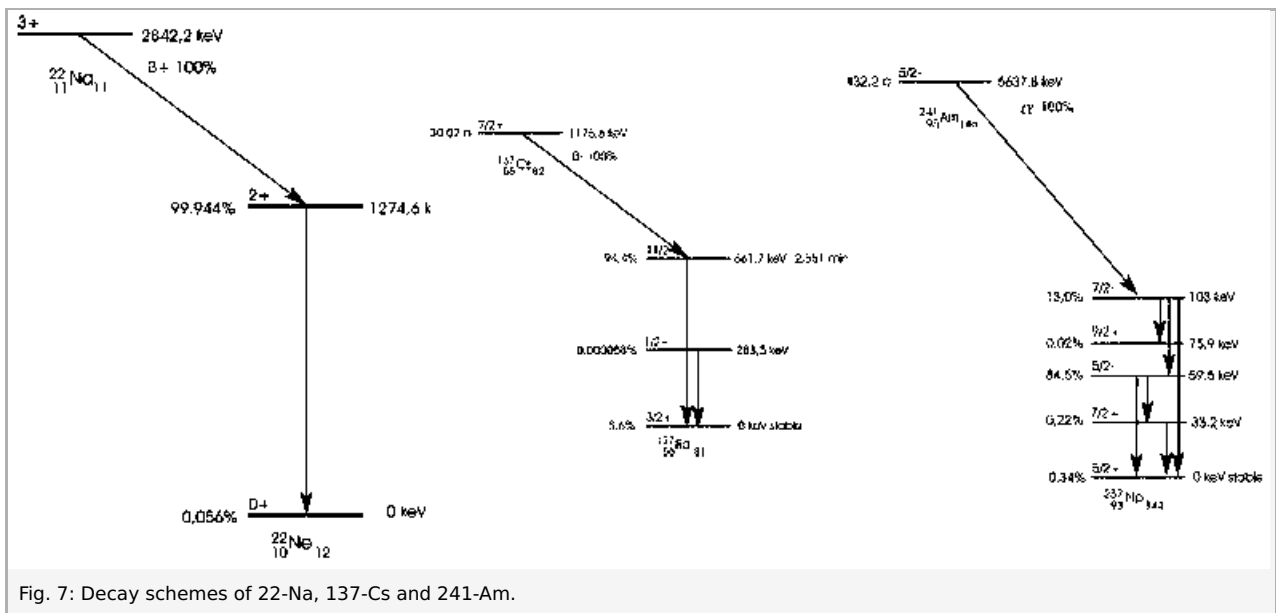


Fig. 7: Decay schemes of 22-Na, 137-Cs and 241-Am.

For each of the processes the fraction of γ -quanta undergoing an interaction per passed layer thickness Δx is constant, it's a reaction probability μ per layer thickness Δx for a fixed γ -energy. For small Δx the absorbed part I_a of the initial intensity (proportional to the number of γ -quanta per unit time) is

$$I_a = \mu \cdot I(x) \cdot \Delta x$$

and

$$I(x + \Delta x) = I(x) - I_a = I(x) - \mu \cdot I(x) \cdot \Delta x$$

thus

$$\frac{dI}{dx} \lim_{\Delta x \rightarrow 0} \Delta x = \frac{I(x + \Delta x) - I(x)}{\Delta x} = -\mu \cdot I(x)$$

is the differential equation describing the absorption. With the border conditions $I(\infty) = 0$ (after an infinitely thick layer there is no intensity left) and $I(0) = I_0$ this has the solution

$$I(x) = I_0 \cdot e^{-\mu \cdot x}$$

The total absorption coefficient $\mu(E)$ is seen as sum of absorption coefficients $\mu_{ph}(E)$ of the photoelectric effect, $\mu_{co}(E)$ of the Compton effect and $\mu_{pa}(E)$ of the pair production effect on Fig. 6. In this experiment, the total absorption coefficient is measured as the slope of a semi logarithmic plot of the intensity vs. the absorption layer thickness for a fixed energy.

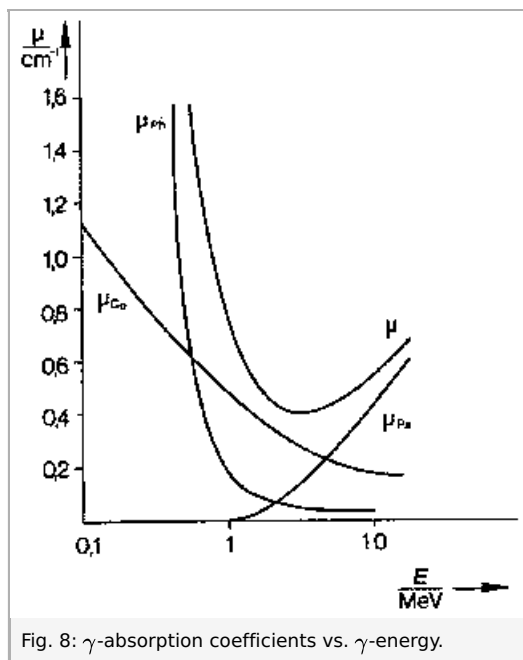


Fig. 8: γ -absorption coefficients vs. γ -energy.

Use the "Regression" tool to determine the base of the exponential function in the recorded measurements. The number of shown decimal digits can be altered with the "Information" button. In Fig. 7 the base reads 0.619646 (± 0.01) thus from

$$I_0 \cdot e^{-\mu \cdot x} = I_0 \cdot b^x$$

follows

$$\mu = -\ln(b),$$

here

$$\mu = (0.48 \pm 0.015) \text{ cm}^{-1}$$

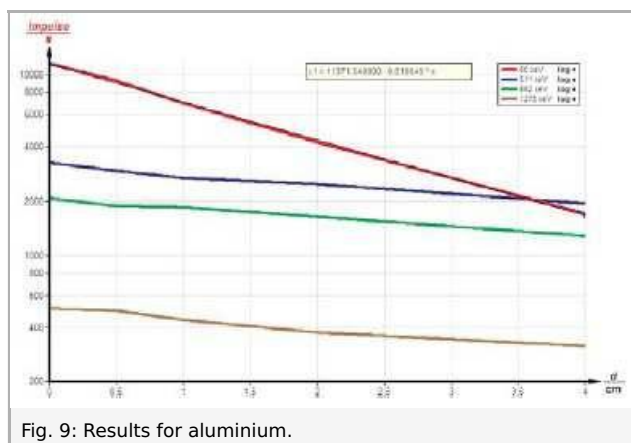


Fig. 9: Results for aluminium.

Plot μ in dependence of the energy. Fig. 8 shows a possible result.

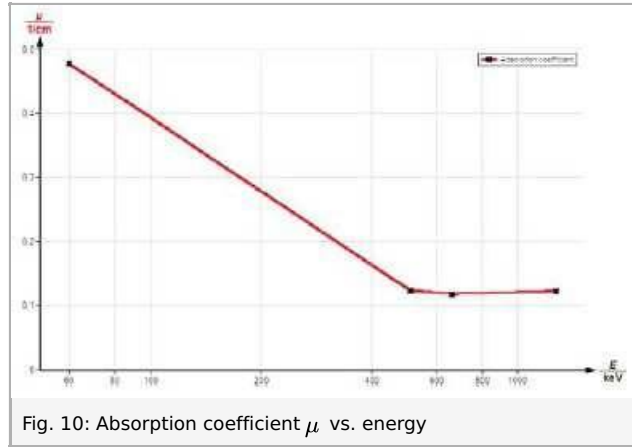


Fig. 10: Absorption coefficient μ vs. energy